



Name: Thymosin beta 4

Pump A: 0.1%Trifluoroacetic in 100% water

Pump B: 0.1%Trifluoroacetic in 100% acetonitrile

Total Flow: 1ml/min

Wavelength: 214nm


Analytical column type: Kromasil 100-5-C18 (4.6*250mm*5 μm)

Dissolution method: 100%H₂O

Inj. Volume: 30ul

Time	Module	Action	Value
0.01	Pumps	B.Conc	15
30.00	Pumps	B.Conc	45
33.00	Pumps	B.Conc	100
38.00	Pumps	B.Conc	100
40.00	Pumps	B.Conc	15
50.00	Controller	Stop	

Chromatogram

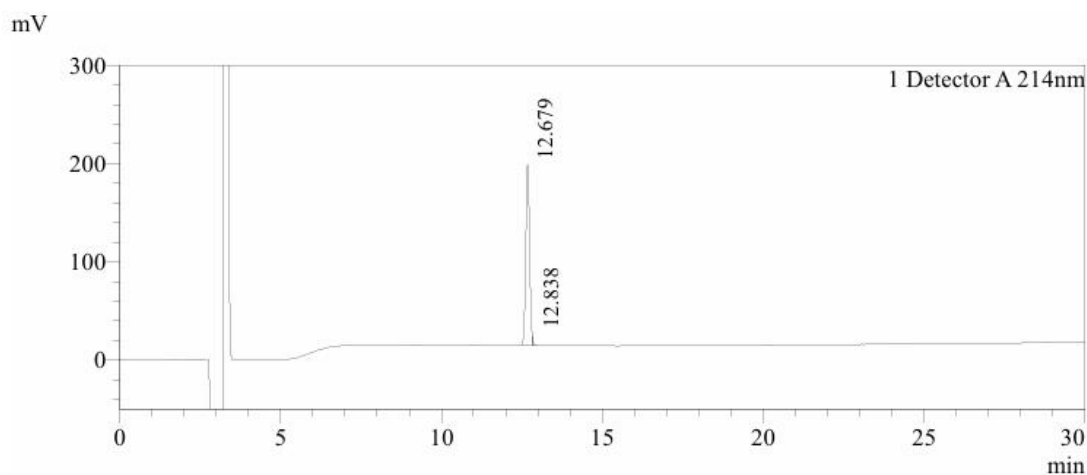

Peak Table

Peak#	Ret. Time	Area	Height	Area%
1	12.679	1521138	183697	99.252
2	12.838	11472	6542	0.748
Total		1532609	190238	100.000

Quality Assurance Officer: 李强

Corrector: 佟丽华

Analyst: 佟琼

